## Ion Cleansing

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
## NTE20

### Table of Contents

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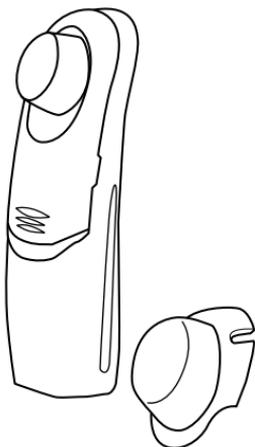Thank you for your purchase. Before using this product, be sure to read the manual and use the product correctly. After reading, keep the manual in a location where users of the product can find it easily.

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Safety Warnings	1~3
Names and Functions of Parts	4
How to Charge	5
Before Using	6・7
How to Use	8~11
When to Charge/ When Using Overseas/ When Using while Plugged In	11
Care and Maintenance	12
Troubleshooting	13・14
Specifications	14

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『Inspecting Ion Cleansing that has been used for a long time!』 ●Make sure the head is not damaged before using.

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# Safety Warnings

- Make sure to read these safety warnings before using, and use the product correctly.
- The warnings listed here are to ensure the product is used properly and safely, and avoid any injury to the user or others or damage to property in advance. Make sure to follow the warnings.
- Safety warnings are separated into categories as follows.

## Danger

For actions that have a significant chance to cause serious injury or death if done improperly.

## Alert

For actions that could cause serious injury or death if done improperly.

## Warning

For actions that could cause injury or property damage if done improperly.

## Examples of Symbols

-  This symbol indicates an action that must not be performed.
  -  (The image to the left indicates "no dismantling.")
-  This symbol indicates an action that must be performed.
-  (Always perform an action indicated by the image to the left.)

## Danger

Never use along with the following medical electronic devices.

- Pacemakers or other medical electronics implanted in the body
- Electrocardiograph or other attachable medical electronics

Prohibited

There is a risk of the medical electronics malfunctioning.

## Alert

Consult with a doctor before use if any of the following apply.

- Heart problems
- Acute (painful) illness
- Fever
- Abnormal blood pressure
- Infectious disease
- Tuberculosis
- Malignant tumor
- Facial neuralgia
- Pregnancy
- While on period
- Atopic eczema
- Allergies
- Particularly sensitive skin
- Hemophilia

Required

Required

- Receiving dental treatment
  - Capillary dilatation due to long-term steroidal hormone treatments or liver dysfunction
  - While under the influence of drugs or alcohol, or extremely tired
- There is a risk of accident or injury.

**Do not use on the following areas.**

- Areas where plastic surgery was performed
  - Eyeballs and upper eyelids
  - Injuries
  - Facial melanosis
  - Inside the mouth or other mucous areas
  - Areas inflamed with pimples or other breakouts
  - Areas with impairment of pain sensation or of perception
  - Areas with metal, plastic, silicone or other objects under the skin
  - Areas that are itchy, have a rash or other stains with physical stimuli
  - Areas inflamed due make-up or other reasons
  - Around the eyes and ears
- There is a risk of accident or injury.

Prohibited

**Do not use on areas other than the face.**

There is a risk of accident or injury.

## Safety Warnings

Do not use the product on wounds or inflammation of skin. It can cause skin problems or pain.

Do not use the product for uses other than those listed in the manual.

It can cause accidents or malfunction, or damage the product.

There is a risk of the power cord disconnecting when used improperly or for long periods, so do not use the product in the following situations:

- If the power cord is twisted and does not return properly.
- If the base of the power cord or adapter has warped or is damaged.

If the power cord is disconnected, there is a risk of sparks occurring and causing burns or electric shock.

Prohibited

Do not damage the power cord. (Do not twist, bend, or pull it unnecessarily, do not place heavy objects on it, do not place it near heat, do not bundle it, do not squeeze it, do not modify it)

There is a risk of fire or electric shock.

Do not use when the power cord is damaged or loose in the outlet. There is a risk of fire or electric shock.

Do not wrap the power cord around the power adapter.

There is a risk that the power cord will disconnect and cause fire, electric shocks or burns.

Do not let children use the product. Do not use in reach of toddlers.

There is a risk of accident or injury.

Do not place it in reach of toddlers. Small parts such as the cap may be swallowed accidentally.

Prohibited

Do not use on toddlers, children or pets.

There is a risk of accident or injury.

Do not use with power sources other than AC100-240V.

There is a risk of fire or electric shock.

If you experience pain or discomfort while using the product or issues with your skin after using the product, stop using the product and consult a medical specialist.

It could cause skin problems.

Required

If the product makes abnormal noises or some other irregularity occurs, stop use immediately.

There is a risk of accident or injury.

Follow directions in manual for frequency and duration of use.

There is a risk of accident or injury.

When first using the product, track condition of skin.


There is a risk of skin problems or poor health.

Remove the power adapter from the outlet when not in use.

There is a risk of fire or electric shock.

Remove  
Power  
Adapter

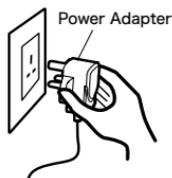No Wet  
Hands

Do not plug in or remove the power adapter with wet hands.

There is a risk of malfunction, electric shock or short circuit.

No Getting  
Wet

Do not wash with water. Do not place the product in water or splash it with water.

There is a risk of malfunction or short circuit.

## Safety Warnings

No Wet  
Areas

Do not use or store the product in a place with high moisture or where it might get wet, such as a bathroom.

There is a risk of malfunction or short circuit.

No  
Dismantling

Do not modify the product. Do not dismantle or repair the product if you are not a repair technician.

There is a risk of fire, electric shock or injury.

### Warning

Do not use at the same time as other beauty goods.

There is a risk of accident or injury.

Do not drop or hit the product. Do not apply strong impacts to the product.

There is a risk of electric shock, fire or injury.

Do not wipe the product with benzene, paint thinner or alcohol, nor spray disinfectant on the product.

There is a risk of cracking, electric shock or fire.

Prohibited

Do not use power adapters other than the one that comes with the product. Do not charge other products with the power adapter.

There is a risk of electric shock or fire.

Do not use for commercial purposes. Do not put it in a place unknown numbers of people may use it.

Continued use can cause malfunctions.

Do not use if damaged.

There is a risk of electric shock, fire or injury.

Prohibited

Do not share with family or strangers.

There is a risk of infection.

Do not place the head directly on skin.

There is a risk of skin problems.

Remove contact lenses and accessories such as earrings and necklaces before using.

There is a risk of accident or injury.

When skin is very rough or the surface is uneven, or when used on areas with hair or whiskers, the product may seem too forceful. In that case, avoid using it on that area.

There is a risk of accident or injury.

Required

Move the cotton around without stopping in one area.

There is a risk of skin problems or health problems.




When removing the power adapter from the outlet, grasp the adapter itself instead of the power cord.

There is a risk of fire or electric shock.

# Names and Functions of Parts

## What is Ion Cleansing?

Tough to clean dirt at the back of pores and old skin cells with negative charges are drawn to the positively charged head and pushed out by electric current. (Note) Not a makeup remover.

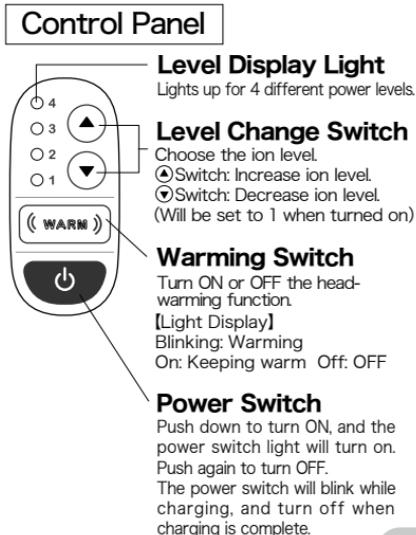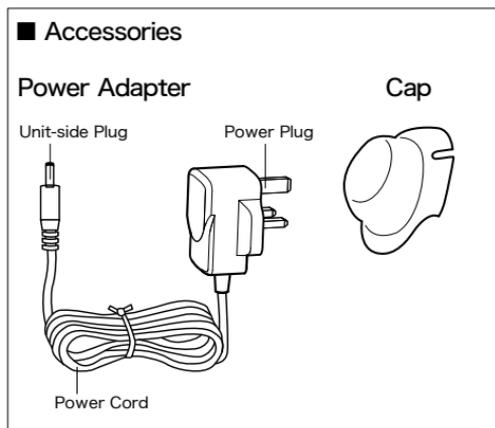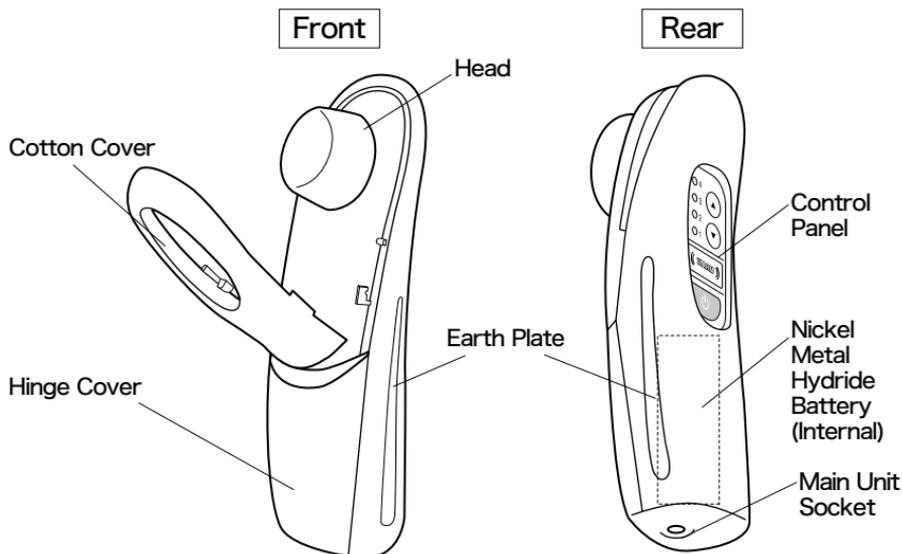\*Follow regional laws when disposing of this product.


# How to Charge

## 1 Confirm the power switch is OFF and connect the power adapter to the main unit and outlet.

- Pushing the unit-side plug of the power adapter too hard or forcefully bending it can damage the main unit socket.

## 2 When the power switch light turns off, remove the power adapter from the main unit and outlet.

- The average charging time is approximately 5 hours.
- While charging, the main unit and power adapter may become warm, but this is not a malfunction.
- One charge allows for using 3 minutes once per day for approximately a month. (Based on factors such as conditions of use and maintenance, the number of uses may decrease.)

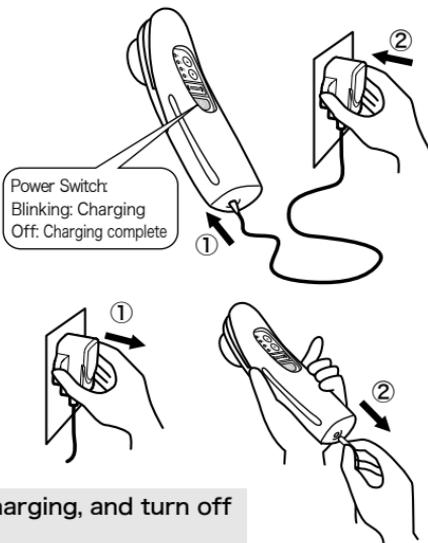The power switch will blink while charging, and turn off when charging is complete.

## How to Have a Long-lasting Nickel Metal Hydride Battery

- When storing the product, fully charge it before storage.
- When not using for a long time, still charge periodically every six months.

When the operating time after charging becomes extremely short, the following procedure can refresh the battery and improve operating time.

1. Push both the power switch and  switch simultaneously for over 10 seconds.
  - When the power switch light, warming switch light and level display lights (1/2/3/4) all start blinking, stop pushing the switches.
  - Let the charge run out until the lights turn off.
2. Charge as normal.

## Warning

- Once charging is complete, remove the power adapter from the outlet.
- Use with an ambient temperature between 10 and 40°C.

# Before Using

## Danger

**Never use along with the following medical electronic devices.**

There is a risk of the medical electronics malfunctioning.

- Pacemakers or other medical electronics implanted in the body
- Electrocardiograph or other attachable medical electronics

## Alert

**Consult with a doctor before use if any of the following apply.**

There is a risk of accident or injury.

- Heart problems
- Acute (painful) illness
- Fever
- Abnormal blood pressure
- Infectious disease
- Tuberculosis
- Malignant tumor
- Facial neuralgia
- Pregnancy
- While on period
- Atopic eczema
- Allergies
- Particularly sensitive skin
- Hemophilia
- Receiving dental treatment
- Capillary dilatation due to long-term steroidal hormone treatments or liver dysfunction
- While under the influence of drugs or alcohol, or extremely tired

**Do not use on the following areas.**

There is a risk of accident or injury.

- Areas where plastic surgery was performed
- Eyeballs and upper eyelids
- Injuries
- Facial melanosis
- Inside the mouth or other mucous areas
- Areas inflamed with pimples or other breakouts
- Areas with impairment of pain sensation or of perception
- Areas with metal, plastic, silicone or other objects under the skin
- Areas that are itchy, have a rash or other stains with physical stimuli
- Areas inflamed due make-up or other reasons

**If you experience any abnormalities with your skin or body while using, immediately stop use and consult with a medical specialist.**

There is a risk of accident or poor health.

## Warning

This appliance is not designed to be used by people (including children) with reduce physical, sensorial or mental capacities, or by people with no prior experience or knowledge, except where they have been able to benefit from supervision or previous instruction relating to use of the appliance from a person responsible for their safety.

## Warning

**Remove contact lenses and accessories such as earrings and necklaces before using.**

There is a risk of accident or injury.

\* Guidelines for skin care are to use once per day, for no more than 3 minutes each time. Using continuously for over 3 minutes can lead to skin problems, so please refrain from doing so.

## 1 Wash your face and clean off any dirt.


- Use cleanser to remove makeup.

Use on skin without any makeup.

## 2 Attach makeup cotton to the main unit.


- ① Open the cotton cover (①, ②) and remove the cap (③).

The cotton cover opens up to 60 degrees. If opened further, the cover will be removed.

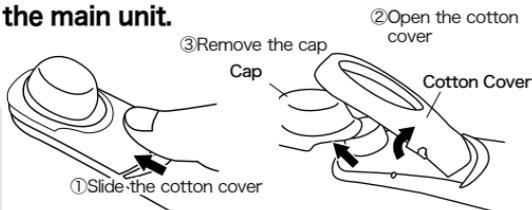
- ② Put cotton on the head (④), and close the cotton cover (⑤).

- Make sure the cotton does not slide off the head.

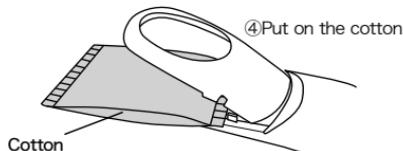- ③ Slide the cotton cover (⑥) and fix the cotton in place.


- Slide the cotton cover down until it hits the hinge cover and there is a clicking noise.

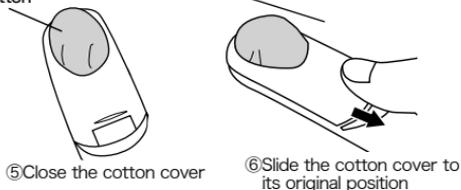Use new makeup cotton each time. Attach it without dipping in beauty wash. (If you dip the cotton in beauty wash before attaching it, then the liquid will spill out when fixing the cotton cover in place)

## Warning

- Do not apply the head directly to the skin. (Always attach cotton.) There is a risk of skin problems.
- Use commercially available makeup cotton. Cotton smaller than the size to the right cannot be attached.
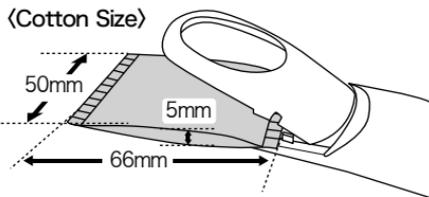
# How to Use

 <b>Alert</b>	<b>Do not use around eyes or ears.</b> (There is a risk of accident or poor health.)
	<b>Do not use on areas where the skin is damaged or injured.</b> (There is a risk of skin problems or pain.)

Use for up to 3 minutes, once per day.  
The power will turn OFF automatically 3 minutes after pushing the power switch.

## 1 Apply beauty wash to the cotton.

Use the beauty wash you normally use or one with cuticle removing effects. (Do not use milky or gel washes)

- Apply enough for the rear of the cotton to absorb it. (Enough that the head can be seen slightly)
- When you apply so much beauty wash that it pours onto the main unit, remove the cotton and wipe off the main unit. (Beauty wash inside the main unit can cause malfunctions.)

Depending on the beauty wash used, effects may be limited.


## 2 Push the power switch on the control panel on the rear of the main unit.

- There will be a beep, display light 1 will turn on, and the warming switch will start blinking and warm the head.
- Approximately 30 seconds later the warming switch will stop blinking and change to a steady light, indicating it is ready. Ion cleansing can be performed even before warming is complete.

## 3 Choose the ion level using the level change switch on the control panel.

- Four ion levels can be chosen, from 1 to 4. The level display light for the level chosen will light up. Choose a level appropriate for your skin (one that does not feel prickly).

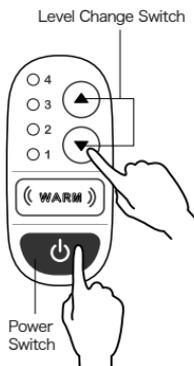
When first using, for people with sensitive skin or people who feel it is strange, it is best to set the level to 1.

## 4 Hold the main unit so your hand is touching the earth plate.

- If your hand is not touching the earth plate during ion cleansing, it will not work.

Make sure that your hand is touching either the left or right earth plate when holding the main unit.

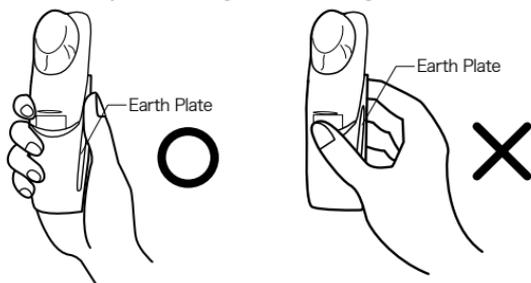
## 5 When the warming switch light turns steady, place the cotton on your skin and slowly move it around.

- Use for up to 3 minutes. The power will turn off automatically 3 minutes after pushing the power switch.

**Warning**  **Move the cotton around, without keeping it on the same spot.**  
(There is a risk of skin problems or poor health.)

- Lightly touch it to skin, and slide it around from the inside to the outside (like the arrows on the image).
- When the cotton loses its moisture or is difficult to slide, apply more beauty wash.
- When you feel minor electric shocks (prickles) while using, apply more beauty wash to the cotton.
- If you still feel prickles after applying more beauty wash, use the level change switch (▼) to decrease the ion level.
- \* If you still feel prickles after taking the above actions, stop using the product.

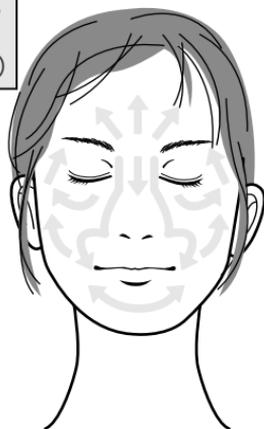
### Warning

- Do not touch the rear control panel during use.  
The power may be turned off or the ion level changes.
- Do not push the control panel switches using nails. This can cause damage to the control panel.
- Use in an ambient temperature between 10 and 40°C. (Suggested ambient temperature when using warming function: 15 to 35°C.)

### 6 It is possible to stop warming the head.

- After turning on the power, the warming switch is set to ON. Warming the skin helps make it easier to remove dirt.
- When warming is not necessary, push the warming switch to stop warming. Pushing again will cause warming to start again.
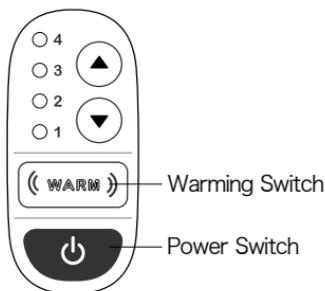
### 7 Push the power switch to turn OFF after using.


- The power will turn off automatically 3 minutes after pushing the power switch. At that time, a beeping noise will alert you.

#### To stop while using

- Push the power switch and the power switch light will turn off.

### 8 Open the cotton cover and remove the cotton.

- Remove the used cotton and wipe off any beauty wash on the main unit with a tissue paper.
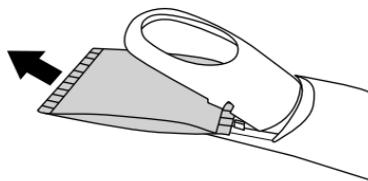
### 9 Put the cap on the head after use and store the product.

- After making sure the head does not become dirty or have moisture, cover the cap.

## Warning

- Use once per day for 3 minutes, and do not use continuously. This can cause skin damage and skin problems.
- After using cotton, dispose of it and do not use again. This can cause skin problems.
- Always wipe off any beauty wash on the main unit. Beauty wash inside the main unit can cause malfunctions.
- After ion cleansing, do normal skin care with beauty wash, lotion, etc.


### Long-term Storage

- Perform maintenance (Page 12) before storing.
1. Put the cap on the head.
  2. Make sure to charge it.
  3. Store in an area free of dust and moisture.

# How to Use

## [When to Charge]

- The signal for low batteries is when the power switch light begins blinking while the power is on and a warning beeping occurs for 10 seconds before the power automatically turns off.  
→ Charge immediately.  
To continue using, use it while plugged in.

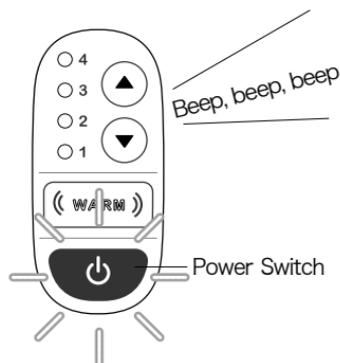
## [When Using Overseas]

- This product can be used with AC100-240V power.


## [When Using while Plugged In]

**1** Connect the power adapter to the main unit and outlet.

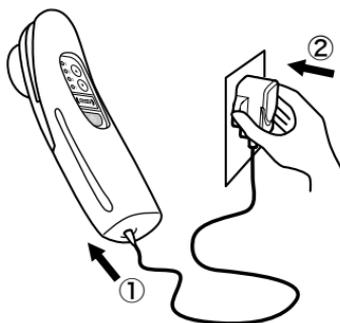**2** Begin cleansing.  
Follow directions on (Pages 8-10) for use.

**3** After use, remove the power adapter from the outlet and main unit.

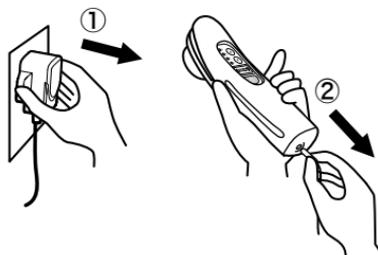**For charging,**  
follow the directions on How to Charge (Page 5).

# Care and Maintenance

Make sure the power is OFF and the power adapter is not attached before performing maintenance.




1 Slide the hinge cover in the direction indicated in the image.

2 Remove the hinge cover (1) and the cotton cover (2).

3 Use a dry, soft cloth to wipe off dirt on the main unit, cotton cover and hinge cover.

● **For difficult to clean dirt**



Soak a dry cloth in soapy water or a neutral cleaning agent cut with water, and wring it out thoroughly before wiping with it.

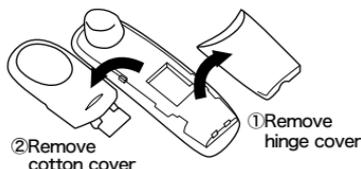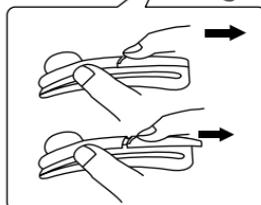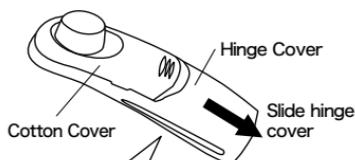4 After attaching the cotton cover to the main unit (1), attach the hinge cover (2).

● Line up the covers with grooves and push them in until there is a clicking noise.

\* Please perform maintenance regularly.
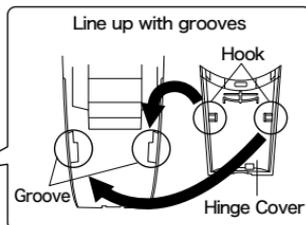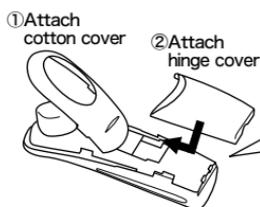
## Warning

- Make sure the power is OFF. If you accidentally turn on the power, it could shock you or cause a malfunction or short circuit.
- Do not wash the main unit with water. Do not immerse it in water or splash water on it. It could cause malfunctions.
- Do not wipe with benzene, paint thinner or alcohol. It could cause deformation or discoloration.

# Troubleshooting

Issues	Possible Reasons	Solutions
Dirt is not cleansing easily.	Your hand is not touching the main unit earth plate.	Hold so that your hand is touching the earth plate.
	The ion level is set too low.	Use the level change switch (⊕) to increase the ion level to a level appropriate for your skin.
	There is not enough moisture in the cotton.	Apply enough beauty wash till you can faintly see the head.
Dirt is hard to see.	There is no visible dirt on your skin.	Continue using as is.
The control panel light is not turning on.	The main unit has no charge.	Charge it or use it while plugged in.
The head is not warming up.	The power switch is not ON.	Turn on the power switch.
	The warming switch is not ON.	Turn on the warming switch.
	Malfunctioning.	Stop use immediately and contact the information number on the product box.
The head is too hot.	The safety device is not functioning properly.	Stop use immediately and contact the information number on the product box.
The power switch begins blinking too quickly.	The product was not used for 3 months or more. (Natural discharge)	Charge the product.
It can only be used a few times after charging.	The product is charged despite there still being charge remaining. (Top-up charge)	Use the battery fully before charging.
	The nickel metal hydride battery is getting old.	Contact the information number on the product box.
Charging does not work.	The power adapter is not connected properly.	Make sure the power adapter is securely connected to the main unit and the outlet.
	The power adapter is malfunctioning.	Contact the information number on the product box.

# Troubleshooting

Issues	Possible Reasons	Solutions
The power cord or power adapter is abnormally hot.	The power adapter is not connected properly.	Firmly connect the unit-side plug of the power adapter to the main unit socket.
	The plug is not connected to the outlet firmly.	The issue may be with the outlet. Contact a local electronics store.
The product cannot be used while the power adapter is plugged in. (It does not turn on even while plugged in.)	The power adapter or main unit is malfunctioning.	Contact the information number on the product box.
When charging or using while plugged in, there is a shock when touching the main unit.	This is due to static electricity build-up in your body.	This is not a malfunction. No action needs to be taken.

In the case of other errors or malfunctions, immediately stop use and contact the information number on the product box.

Specifications	Product Name	Ion Cleansing	Rated Time	3 minutes (auto OFF)	
	Product Number	NTE20	Charge Time	Approx. 5 hours	
	Power Source	Power Adapter	AC100-240V 50/60Hz	Weight	110g (main unit)
		Main Unit	Nickel Metal Hydride Battery (DC3.6V, 750mAh) 1 unit		
	Power Consumption	When plugged in 10W		Accessories	Cap Power Adapter
Dimensions	Height 130mm x width 41mm x depth 38mm (main unit)				

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